PRINCESS®

Gebruiksaanwijzing
Instructions for use
Mode d'emploi
Anleitung

TYPE 541000
LICE COMB

NL	Gebraiksaanwijzing	4
UK	Instructions for use	5
F	Mode d'emploi	6
D	Anleitung	8

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PRINCESS LICE COMB

ART. 541000




Read these instructions before use and keep them for future reference.

The PRINCESS LICE COMB is a chemicals-free lice treatment device which uses the latest technology.

The dual-combs move in opposite direction from each other and automatically remove most lice during treatment. For a quick solution to the problem, it is recommended to use the comb daily in order to remove both lice and their eggs (nits).

BEFORE USE

- Remove all packaging.
- Use a coin to open the battery compartment cover by turning it anti-clockwise (fig 1).
- Insert 4 x AA 1.5V Alkaline batteries as shown in fig 2.
- Use a coin to close the battery compartment cover by turning it clockwise (fig 3).

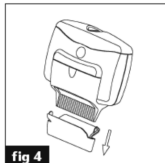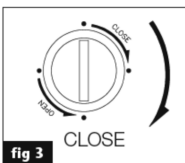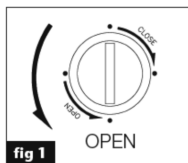USE

- Place a towel or sheet over the shoulders.
- Remove the transparent protective cover (fig 4).
- Push the '0/1' slide to the '1' position (fig 5).
- Start combing the hair with the appliance, working from the root to the tip of the hair. Pay extra attention to those areas where lice are most likely to occur i.e. near the scalp around the ears and at the back of the neck.

For best results, first wash the hair using shampoo, then condition with a conditioner, and then comb or brush it. Finally comb with the PRINCESS LICE COMB while the hair is still damp.

Combing should be done in slow strokes while holding the comb at an angle to the scalp (fig 6).

Rinse off any lice that stick to the comb at the end of a stroke.


Each area should be combed at least 3 times.

- As nits hatch within 7 to 9 days, daily treatment should continue for 12 to 14 days to ensure that all lice are removed.
- While combing long hair, the comb may get stuck. In that case, remove the comb from the hair and resume combing about 1" (2.5mm) back.
- At the end of the treatment, push the '0/1' slide to the '0' position to switch off the appliance (fig 7).
- Clean the unit (see "CLEANING" section).
- Replace the protective cover on to the main unit.

Note:

It is recommended to use a separate comb for each person being treated.

Remove the batteries from the main unit if you are not going to use the comb for a long period of time.

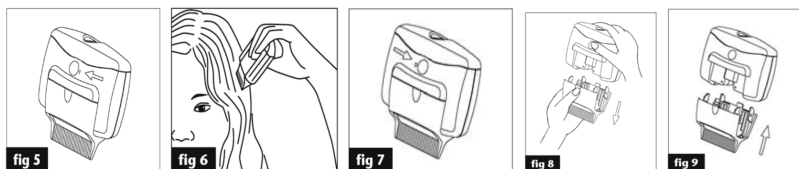CLEANING

Turn off the appliance before cleaning.

Remove the comb assembly from the main unit (fig 8).

Wash the comb assembly in hot soapy water (65°C/150°F) and then rinse.

Use a soft damp cloth to clean the main unit.

Dry the comb assembly completely and then replace it on to the main unit (fig 9).

Put the protective cover back on and store.

PRINCESS TIPS

- This appliance is meant for household use only. Do not use the appliance outdoors.
- The appliance can be used on wet or dry hair.
- Do not use your lice comb on hair that has been treated with gel or any other hair styling product.
- This comb is intended for use on scalp hair only.
- Avoid combing near the eyes.
- Do not use in the shower.
- Do not immerse the main unit in water or any other liquid.
- Children should be supervised while using the appliance.
- The use of this appliance by children or adults with a physical, sensory, mental or motorial disability, or lack of experience and knowledge can be hazardous. Persons responsible for their safety should give explicit instructions and supervise their use of the appliance.

BATTERIES

Batteries constitute chemical waste and may not be discarded in the normal fashion. Put empty batteries in the recycling bin or take them separately to your local refuse department following local procedures.

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